

\$

Save Time

Cost Effective



Easy-to-use

TECHNICAL SPECIFICATIONS

 System Compatibility:
 NMReady-60e or 60PRO

 Dimensions:
 10.3x10.3x14.2"/26x26x36cm

 Weight:
 12 lbs/5.5 kg

 Voltage:
 24 V

 Time per tube change:
 ~10sec

 Temp:
 32°C

 Carousels:
 2

CONTACT US



Nanalysis Corp 1-4600 5th St. NE Calgary, AB, Canada T2E 7C3 Nanalysis GmbH Haydnplatz 3 D-76133 Karlsruhe, Germany

49.159.04529975

1.855.NMReady




sales@nanalysis.com



nanalysis



AUTOsample-60

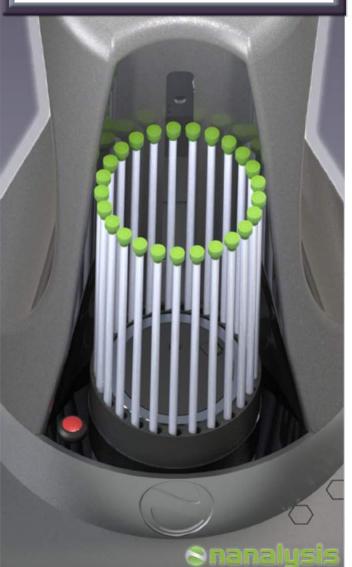

Snanalysis

Maximize your efficiency by automating routine work.



Introducing the affordable, easy-to-use AUTOsample-60.

The AUTOsample-60 mounts on top of any NMReady-60e or 60PRO for painless interconversion between autosampler and manual mode, while still maintaining the notably minimal footprint. Data acquisition and processing is controlled at the touch of a finger through the onboard software.

EASY CONFIGURATION

Simple mounting on top of NMReady-60

- Operated through onboard
 software via NMReady USB connectors
- Interchangeable 25 NMR tube
 carousels (standard, constricted or J-Young tubes)
 - Robust and reliable tube changes ensure safety and accuracy
 - Optional sample tempering

Emergency stop

Interchangeable 25 tube capacity carousel

SIMPLE OPERATION

Use onboard software to configure data acquisition of up to 25 different samples at one time.

Set up structure elucidation or conformation protocols per tube e.g., ¹H, COSY, HSQC, ¹³C{¹H}
 QA/QC for series of

samples

Temperature ---- controlled samples

Sturdy, no-slip base

- Configurable experiment protocols and tube configuration

Easy connection to onboard software